

# CA100 Wall Mounted – Right Hand Sink CB100 Wall Mounted – Left Hand Sink CC100 Wall Mounted – Center Sink CC300 Wall Mounted – Dual Approach

## Autopsy Sink



CA100 & CB100 with DC100 Cart





CC100



CC300

*Photos may include optional features that are not standard*

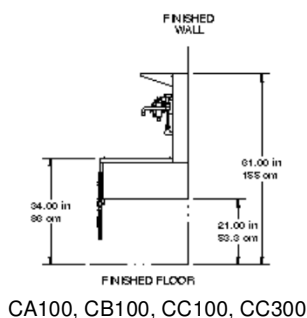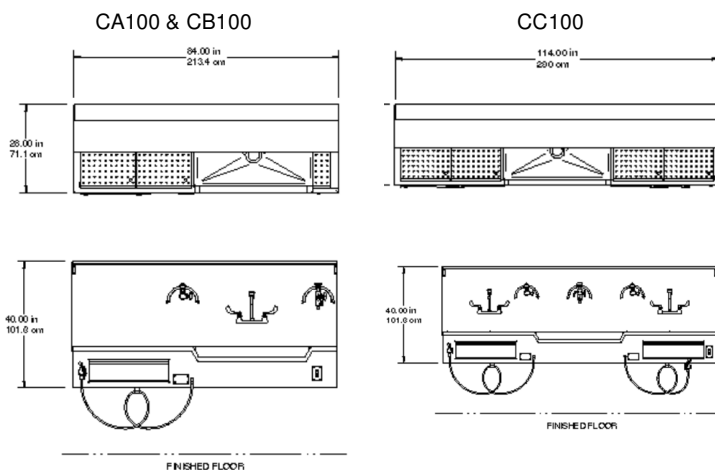

CA100, CB100, CC100, CC300

### Standard Features

- 84"L x 28"W x 61"H – **CA100 & CB100**
- 114"L x 28"W x 61"H – **CC100**
- 150"L x 28"W x 61"H – **CC300**
- Work surface 34" from finish floor to bottom of work surface
- Hydro aspirator with reversing flow valve
- Hand Spray Assembly with durable 8' high pressure hose.
- Single lever hot/cold water mixing valve with write blade handles
- Interchangeable Grid Plates
- Rinse basket grid plate.
- GFCI duplex outlets with waterproof covers
- Instrument drawer flush to the table
- Autopsy Cart latch – compatible with DA000 & DB000 cart chassis only

### Optional Features

#### Disposal Options

BL800 Disposal 1 HP

#### Exhaust Options

CO007 Exhaust Duct

#### Rinse Options


CO003 Sink Perimeter Rinse Assembly  
CO020 Organ Rinse Basket  
CO022 Evidence Screen

#### Attachment Accessory Options

CO001 Autopsy Sink Leg Frame – **CA100 & CB100 Only**  
CO002 Autopsy Sink Leg Frame – **CC100 Only**  
CO028 Autopsy Sink Leg Frame – **CC300 Only**  
CO004 Gas Valve Splash Mounted  
CO005 Air Valve Splash Mounted  
CO006 Safety Eye Wash Splash Mount  
CO008 LED Light Assembly  
CO012 Inch/Centimeter Ruler  
CO016 Scale stand, wall mounted  
CO017 Scale stand, table mounted  
SA100 Seismic Anchor Kit (*For use with CO001*)

#### Assembly is required.

Mopec installation is available upon request.

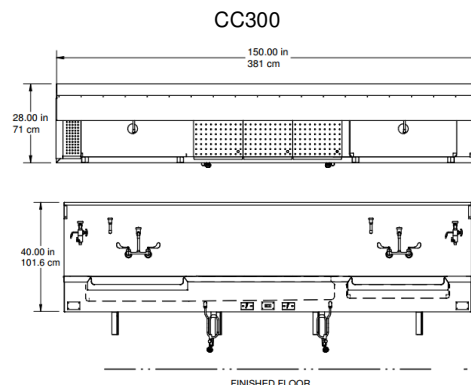

# CA200 Wall Mounted Elevating– Right Hand Sink CB200 Wall Mounted Elevating – Left Hand Sink CC200 Wall Mounted Elevating – Center Sink

## Autopsy Sink





CA200 & CB200 with DC100 Cart




CC200

*Photos may include optional features that are not standard*

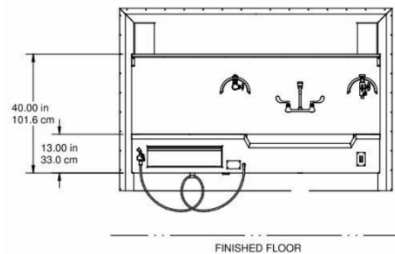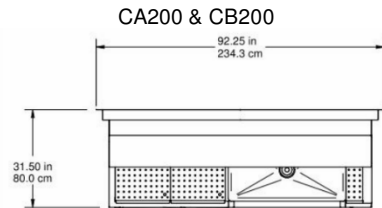

CC200

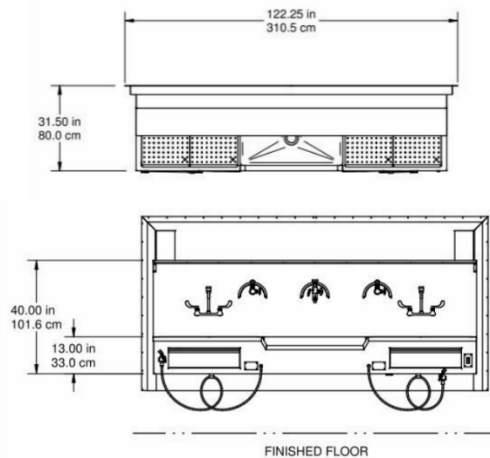

### Standard Features

- 84"L x 28"W x 58.5"H to 66.5" H – **CA200 & CB200 Only**
- 114"L x 31.5"W x 58.5"H – **CC200 Only**
- Stainless Steel Sink
- Hydro aspirator with reversing flow valve
- Hand Spray Assembly with durable 8' high pressure hose.
- Single lever hot/cold water mixing valve with wrist blade handles
- Interchangeable Grid Plates
- Rinse basket grid plate.
- GFCI duplex outlets with waterproof covers
- Instrument drawer flush to the table
- Autopsy Cart latch

### Optional Features

#### Disposal Options

BL800 Disposal 1 HP

#### Exhaust Options

CO018 Ventilation Grill w/Exhaust Duct

#### Rinse Options


CO003 Sink Perimeter Rinse Assembly  
CO020 Organ Rinse Basket  
CO022 Evidence Screen

#### Attachment Accessory Options

CO004 Gas Valve Splash Mounted  
CO005 Air Valve Splash Mounted  
CO006 Safety Eye Wash Splash Mount  
CO008 LED Light Assembly  
CO012 Inch/Centimeter Ruler  
CO016 Scale stand, wall mounted  
CO017 Scale stand, table mounted

**Assembly is required.**

**Mopec installation is available upon request.**

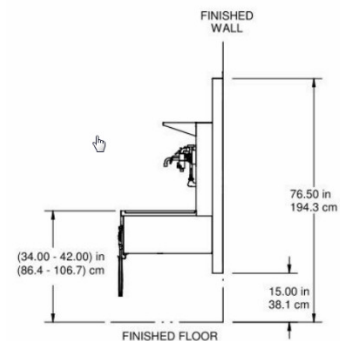



# CE800 Autopsy Sink Service Island




## Autopsy Sink



CE800 Autopsy Service Island  
shown with DC100 Autopsy Cart



*Photos may include optional features that are not standard*

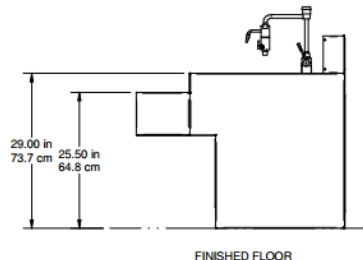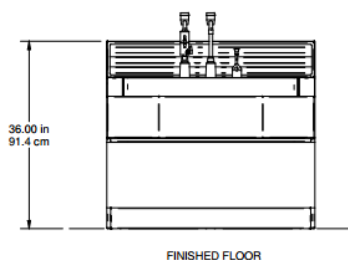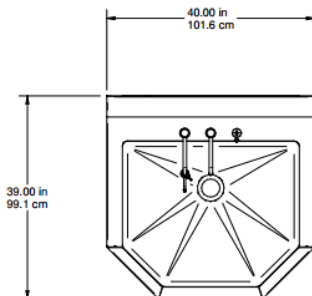

### Standard Features

- 40"L x 39"W x 36"H
- Stainless Steel
- Aspirator collection: device consists of two (2) 4 liter heavy duty polypropylene bottles constructed for heavy vacuum application.
- Bottles are autoclavable and chemical resistant.
- Stainless steel carriage for bottles with quick disconnect racket designed to be mounted anywhere on the workstation.
- Aspirator hose provided with hose caddy.
- Swing spout faucet provides single lever handle hot/cold mixing valve.
- Deck mounted fixtures capable of withstanding formalin and saline environment.
- Hydro aspirator located on table near sink basin.
- Workstation ventilated with a wash down ventilation grill. Ventilation is routed through the base of the unit and into the in-house ventilation system.

### Optional Features

#### **Disposal Options**

BL800 Disposal 1 HP

#### **Rinse Options**

CO003 Sink Perimeter Rinse Assembly

CO009 Perimeter Rinse

CO036 Distilled Water Fixture

CO037 Hand Spray Rinse

#### **Attachment Accessory Options**

CO017 Scale stand, table mounted

CO012 Inch/Centimeter Ruler

BK610 Fluid Collection System

#### **Formalin Options**

CO035 Formalin Fixture

**Assembly is required.**

**Mopec installation is available upon request.**





# CE850 Autopsy Sink Service Island


## Autopsy Sink



*Photos may include optional features that are not standard*

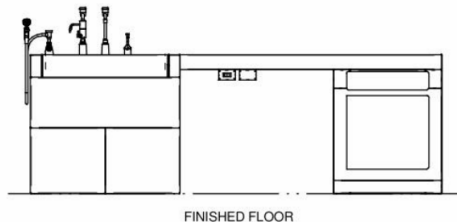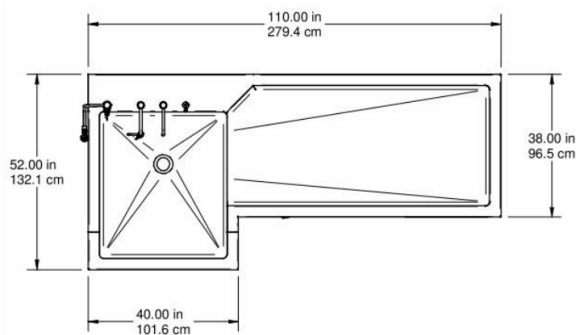

FINISHED FLOOR

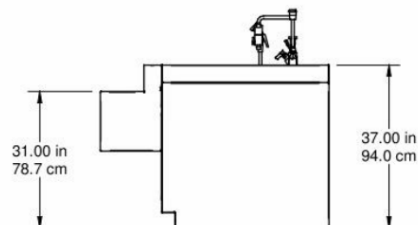

FINISHED FLOOR

### Standard Features

- 100"L x 52"W x 37"H
- Stainless Steel
- Service Table top is 70"L x 38" W
- Stainless Steel Sink
- Hydro aspirator with reversing flow valve
- Hand Spray Assembly with high pressure hose.
- Single lever hot/cold water mixing valve.
- Interchangeable Grid Plates
- Rinse basket grid plate.
- GFCI duplex outlets with waterproof covers
- Instrument drawer flush to the table

### Optional Features

#### **Disposal Options**

BL800 Disposal 1 HP

#### **Rinse Options**

CO023 Table Rinse Control Valve with Nozzles

#### **Attachment Accessory Options**

CO017 Scale stand, table mounted

CO012 Inch/Centimeter Ruler

BA078 Headrest, stainless steel

BA025 Headrest, rubber

BE120 Autopsy Table Head Rinse (**Must use with BA025 and BA078**)

**BK610 Fluid Collection System**

**Assembly is required.**

**Mopec installation is available upon request.**



# CH100 Scrub Sink – Single Faucet

# CH200 Scrub Sink – Double Faucet

---

## Scrub Sink

### Standard Features

#### **CH100**

- 30" (76.2cm)L x 15 3/8" (39cm)W x 23" (58.42cm)H
- Sink dimensions 26 3/4" (68cm)L x 15 3/8" (39cm)W x 11" (28cm)H.

#### **CH200**

- 48" (121.92cm)L x 19" (48.26cm)W x 23" (58.42cm)H
- Sink dimensions 44 3/4" (113.6cm)L x 15 1/3" (39.05cm)W
- Stainless steel scrub sink with 12" backsplash
- Hot & cold foot pedal valve with gooseneck faucet



CH100



CH200





**Assembly is required.**  
**Mopec installation is available upon request.**

*Photos may include optional features that are not standard*

CH100

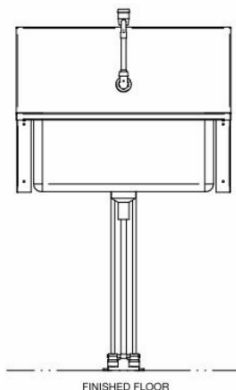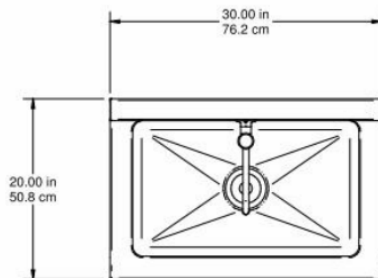

CH200

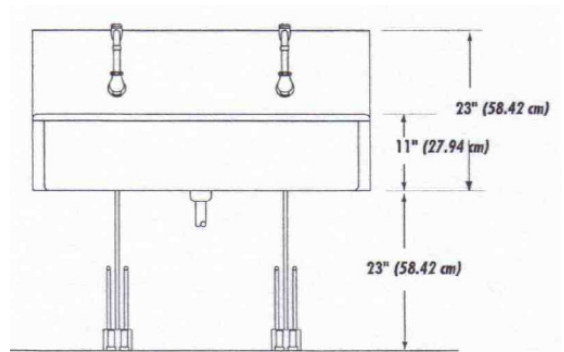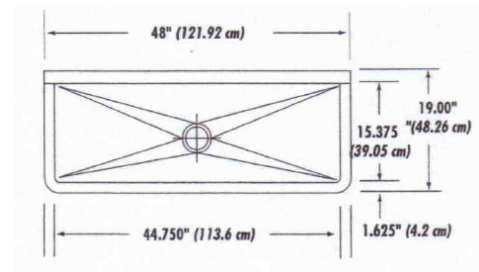




# CF100 – Embalming Sink – Right


# CF200 – Embalming Sink – Left

# CF300 – Embalming Sink - Center

## Embalming Sink




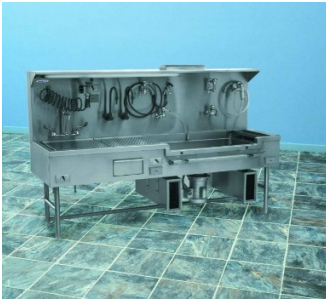


CF100 & CF200

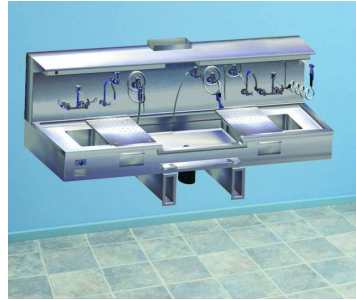

CF300

*Photos may include optional features that are not standard*

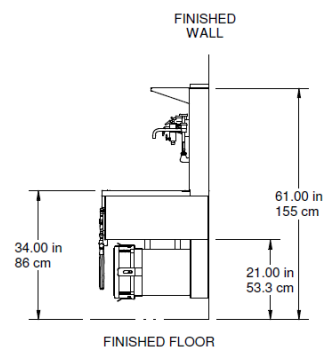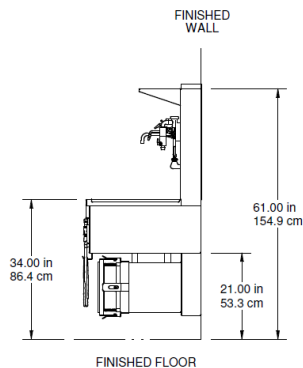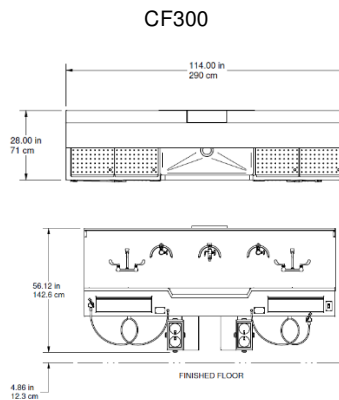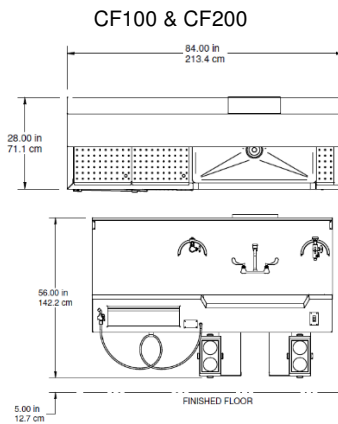

### Standard Features

- 84"L x 28"W x 56"H – **CF100 & CF200**
- 114"L x 28"W x 56.12"H – **CF300**
- Stainless steel embalming sink: 12"L x 12"W x 9"Deep
- Hot/Cold Water Fixture with Wrist Blade Handles
- Gooseneck faucet and vacuum breaker fixture
- Sink spray assembly with high pressure hose – vacuum broke hydro aspirator.
- Removable stainless steel perforated grid plates – Two (2)
- Built-in sponge bowl-stainless steel drawers with drip rail handle.
- Quick Disconnect Air Ventilation Chambers – Removable Top Shelf
- GFCI duplex outlets.
- Body Rinse consists of cold water valve, serrated hose fitting, vacuum breakers and 8" clear tubing.

### Optional Features

#### Disposal Options

BL800 Disposal 1 HP

#### Exhaust Options

CO007 Exhaust Duct

#### Rinse Options

CO003 Sink Perimeter Rinse Assembly

CO020 Organ Rinse Basket

CO022 Evidence Screen

#### Attachment Accessory Options

CO001 Autopsy Sink Leg Frame – **CF100 & CF100 Only**

CO002 Autopsy Sink Leg Frame – **CF300 Only**

CO004 Gas Valve Splash Mounted

CO005 Air Valve Splash Mounted

CO006 Safety Eye Wash Splash Mount

CO008 LED Light Assembly

CO012 Inch/Centimeter Ruler

CO016 Scale stand, wall mounted

CO017 Scale stand, table mounted

CO027 Switched Power Source for Embalming Equipment

**Assembly is required.**

**Mopec installation is available upon request.**

# FD200 – Embalming Sink

---

## Embalming Sink



### Standard Features

- 32"L x 35.3"W x 58"H
- Large sink: 17"L x 25"W x 10"Deep
- Built in sterilization tray
- Deluxe hydro aspirator
- GFCI duplex outlet
- Heavy duty table top with positive 2" marine edge
- Heavy crease to drain

### Optional Features

#### **Disposal Options**

BL800 Disposal 1 HP

#### **Rinse Options**

CO009 Perimeter Rinse Assembly

#### **Attachment Accessory Options**

CO024 Large Fitted Sink Cover for Aerosols Containment



CO024 Optional Fitted Sink Cover shown



FD200 shown with DC100 Autopsy Cart

*Photos may include optional features that are not standard*

**Minor setup is required.**

**Mopec installation is available upon request.**