

## cellIQ™

CO<sub>2</sub> Incubator

8.1 cu.ft. | 230 L



MCO-230AICUVL-PA



Enhance the production of your cell culture experiments with our large capacity CO<sub>2</sub> incubator. Ideal for high throughput culturing or protocols that may require the use of experimental apparatus.

### Precise and Repeatable Environment


The patented Direct Heat and Air Jacket system works with the IR2 sensor to provide precise temperature and CO<sub>2</sub> levels, respectively. InCu-saFe® and SafeCell™ UV function together to prevent contamination.

### Time-Saving Decontamination

The optional high speed decontamination system uses hydrogen peroxide (H<sub>2</sub>O<sub>2</sub>) vapor to safely clean the chamber in less than three hours while providing a kill rate with at least a six log reduction of major contaminants.

### Improved Use and Maintenance

A color LCD touchscreen panel delivers full control, even with gloved hands, while a standard USB port makes log data transfer to PC convenient. Integrated shelf channels make cleaning easy within the large interior and allow more capacity for cell culture containers and apparatus.

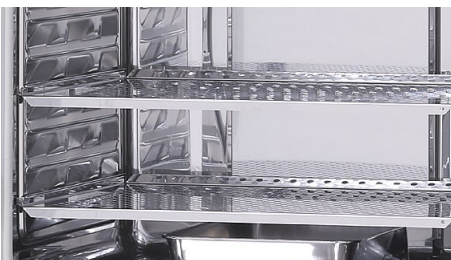

### Germicidal Interior

InCu-saFe copper enriched, stainless steel alloy offers germicidal properties of elemental copper without corrosion or discoloration.



### Contamination Control

The isolated UV lamp decontaminates circulated air and the humidity water reservoir without harming the cultured cells.



### Intuitive Control

The control panel interface allows for easy programming of temperature, CO<sub>2</sub> and other internal conditions of your incubator.




Model MCO-230AICUVL-PA

### Direct Heat and Air Jacket System

Accurate temperature levels are regulated by the Direct Heat and Air Jacket system which functions under microprocessor control.


### Direct Heat and Air Jacket System

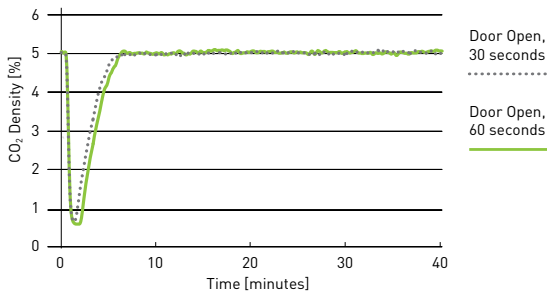

Proprietary IR2 sensor achieves ultrafast CO<sub>2</sub> recovery without overshoot, even following multiple door openings.

### Dimensions

Unit : inch | mm

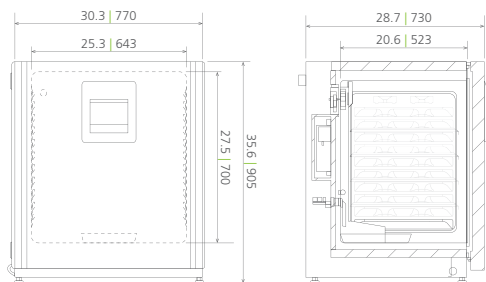

**PHC Corporation of North America**  
 1300 Michael Drive, Suite A, Wood Dale, IL 60191  
 Toll Free USA (800) 858-8442, Fax (630) 238-0074  
[www.phcd.com/us/biomedical](http://www.phcd.com/us/biomedical)

Model Number		MCO-230AICUVL-PA	With Optional H <sub>2</sub> O <sub>2</sub> Decontamination
External Dimensions (W x D x H) <sup>1)</sup>	inches   mm	30.3 x 28.7 x 35.6	770 x 730 x 905
Internal Dimensions (W x D x H)	inches   mm	25.3 x 20.6 x 27.6	643 x 523 x 700
Volume	cu.ft.   liters	8.1   230	
Net Weight	lbs   kg	198   90	
<b>Performance</b>			
Warranty <sup>2)</sup>		3 years parts and labor, 5 years CO <sub>2</sub> sensor	
Temperature Control Range and Fluctuation <sup>3)</sup>	°C	+5 above ambient to +50, ± 0.1	
Temperature Uniformity <sup>3)</sup>	°C	± 0.25	
CO <sub>2</sub> Control Range and Fluctuation <sup>3)</sup>	%	0 to 20, ± 0.15	
Humidity Level & Fluctuation	% RH	95 at 37°C, ± 5	
<b>Control</b>			
Controller		Microprocessor	
Temperature Sensor		Thermistor	
Display		Color LCD touchscreen readable in 0.1 increments	
Sensor	CO <sub>2</sub>	Dual filter IR	
<b>Construction</b>			
Exterior Material		Painted steel (rear cover coated steel)	
Interior Material		Stainless steel copper enriched alloy	
Outer Door	qty	1	
Field Reversible Door		Included	
Inner Door	qty	1 (sealing tempered glass with positive latch)	
Humidity Pan	qty	1 (stainless steel)	
Shelves	qty	4 (stainless steel copper enriched alloy)	
Shelf Dimension (W x D x H)	inches   mm	24.4 x 17.7 x 0.5   620 x 450 x 12	
Max. Load per Shelf	lbs   kg	15.4   7	
Max. Total Load	lbs   kg	61.6   28	
Max. Shelf Capacity	qty	10	
Access Port / Position	qty	1; rear upper left	
Access Port Diameter	inches   mm	1.2   30 (with silicone (non-VOC) stopper)	
Leveling Feet	qty	4	
<b>Decontamination Control</b>			
InCu-safe Chamber, Plenum, Shelves	passive	Included (stainless steel copper enriched alloy)	
SafeCell UV Light System	passive/active	MCO-170UVS-PA (included)	
Hydrogen Peroxide (H <sub>2</sub> O <sub>2</sub> ) Vapor	active	Optional	Included
<b>Alarms (V=Visual Alarm, Buzzer Alarm, R=Remote Alarm)</b>			
Power Failure		R	
Temperature Deviation	high	V-B-R	
Gas Deviation	CO <sub>2</sub>	V-B-R	
Supply Empty	CO <sub>2</sub>	V-B-R	
Door Open		V-B	
<b>Electrical and Noise Level</b>			
Power Supply		115V, 1Ø, 60Hz, NEMA 5-15P requires NEMA 5-15R receptacle	
Noise Level <sup>4)</sup>	dB(A)	25	
<b>Options</b>			
Hydrogen Peroxide (H <sub>2</sub> O <sub>2</sub> ) Vapor Board		MCO-170HB-PA <sup>5)</sup>	Included
Outer Door – Password Access Electronic Lock		MCO-170EL-PW <sup>5)</sup>	Included
H <sub>2</sub> O <sub>2</sub> Vapor Generator		MCO-HP-PW <sup>5)</sup>	
H <sub>2</sub> O <sub>2</sub> Reagent	package of 6	MCO-H2O2-PV	
CO <sub>2</sub> Gas Pressure Regulator	psi	0 – 15; MCO-100L	
Automatic CO <sub>2</sub> Cylinder Changeover System		MCO-21GC-PW	
InCu-safe Shelf		MCO-230ST-PW	
InCu-safe Shelf-Reinforced		MCO-230RT-PW	
Double Stacking Bracket <sup>6)</sup>		MCO-230SB-PW	
Roller Base		MCO-230RB-PW	
4-20mA Analog Output		MCO-420MA-PW	
<b>Optional Communication System</b>			
Wireless, Cloud-based, Automatic Data Management		LabAlert® Monitoring System	

<sup>1)</sup> Exterior dimensions of main cabinet only, excluding handle and other external projections

<sup>2)</sup> Current warranty offered at time of printing and may be subject to change; US and Canada only

<sup>3)</sup> Ambient temperature 23°C, setting 37°C, CO<sub>2</sub> 5%, no load

<sup>4)</sup> Nominal value, background noise 20 dB(A)

<sup>5)</sup> MCO-230AIC series requires MCO-170HB-PA, MCO-170EL-PW and MCO-HP-PW for H<sub>2</sub>O<sub>2</sub> Decontamination

<sup>6)</sup> If stacking two incubators, make sure the double stacking dedicated secure hardware and spacer are used  
 Note: Additional options available.

Specifications are subject to change without notice.

For latest specification information contact PHC Corporation of North America at [info@us.phcd.com](mailto:info@us.phcd.com).  
 Performance data herein is based on independent testing at time of publication.