


The Plus Series Built-In Undercounter Refrigerators feature a microprocessor temperature controller and display, which allows for precise temperature control, verification and recovery versus mechanical style thermostats. Refrigerators feature an interior fan driven, forced-draft air circulation system provides excellent temperature uniformity and recovery after door openings.

Models are designed with a "front-breathing" ventilation system, eliminating the need for clearance around the top and sides of the unit; enabling the user to incorporate these units into cabinets and other applications where free space is limited.



DAI-HC-UCBI-0404G-ADA

DAI-HC-UCBI-0204

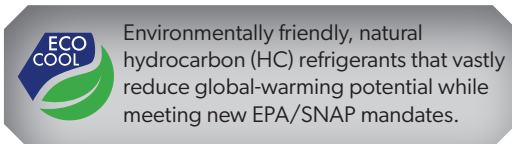





### FEATURING

- Built-in, front breathing technology
- Microprocessor temperature controller
- Adjustable temperature range: 1°C to 10°C
- Cycle defrost
- Remote alarm contacts
- Back wall evaporator cover
- Forced draft circulation
- Keyed door lock
- Magnetic door gasket for positive seal
- Leveling legs
- LED interior lights (select models)
- Probe access port
- ADA compliant units available
- UL/C-UL listed
- Two year parts and labor warranty, plus an additional three year compressor parts warranty



- Digital temperature display
- Audible and visual temperature alarms
- Resettable Min/Max
- °C/°F switchable

Glass Door Model	Ft <sup>3</sup>	Set Point	Defrost	Doors	Shelves	Width"	Depth"	Height"
 DAI-HC-UCBI-0404G*	4.6	4°C	Cycle	1 Glass	2 Adj./1 Fixed	23 ¾	25 ½	33 ¾
 DAI-HC-UCBI-0404G-ADA*	4.6	4°C	Cycle	1 Glass	2 Adj./1 Fixed	23 ¾	25 ½	31 15/16
Solid Door Model								
 DAI-HC-UCBI-0204*	2.5	4°C	Cycle	1 Solid	4 Adj.	17 ¾	22 ¾	30 ¾
 DAI-HC-UCBI-0404SS*	4.6	4°C	Cycle	1 Stainless Solid	2 Adj./1 Fixed	23 ¾	24 ½	33 ¾
 DAI-HC-UCBI-0404*	4.6	4°C	Cycle	1 Solid	2 Adj./1 Fixed	23 ¾	24 ½	33 ¾
 DAI-HC-UCBI-0404-ADA*	4.6	4°C	Cycle	1 Solid	2 Adj./1 Fixed	23 ¾	24 ½	31 15/16

Specifications are subject to change  
\*Left hinge model available