

The Plus Series Freestanding Undercounter/Countertop Freezers feature a microprocessor temperature controller and display, which allows for precise temperature control, verification and recovery versus mechanical style thermostats. **For manual defrost units, periodic manual defrosting is required.**

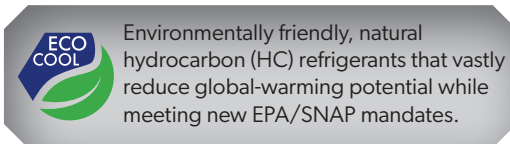


DAI-HC-UCFS-0430






DAI-HC-UCFS-0120A

FEATURING

- Microprocessor temperature controller
- Adjustable temperature range:
(-20°C): -15°C to -25°C
(-30°C): -25°C to -35°C
(-40°C): -35°C to -45°C
- Auto & Manual defrost available
- Remote alarm contacts
- Keyed door lock
- Magnetic door gasket for positive seal
- Leveling legs
- Probe access port
- UL/C-UL listed
- Two year parts and labor warranty, plus an additional three year compressor parts warranty
- These models require 4" of clearance around back and sides of the unit for proper ventilation. Improper installation will void the warranty and will lead to temperature maintenance issues, overheating and possible failure



- Digital temperature display
- Audible and visual temperature alarms
- Resettable Min/Max
- °C/°F switchable

Freezer Model	Ft ³	Set Point	Defrost	Doors	Shelves	Width"	Depth"	Height"
 †DAI-HC-UCFS-0120A*	1.3	-20°C	Auto	1 Solid	2 Adj.	23 ¾	24	20 ¾
 DAI-HC-UCFS-0120*	1.7	-20°C	Manual	1 Solid	2 Adj.	23 ¾	24	20 ¾
 DAI-HC-UCFS-0520	4	-20°C	Manual	1 Solid	2 Fixed + 1 Basket	23 ⅝	24 ¼	32 ¾
 DAI-HC-UCFS-0430	4	-30°C	Manual	1 Solid	3 Adj.	23 ¾	24	32 ⅝
 DAI-HC-UCFS-0440	4	-40°C	Manual	1 Solid	3 Adj.	23 ¾	24	32 ⅝

Specifications are subject to change

*Left hinge model available

†DISCLAIMER: Auto defrost freezers incorporate an electric heater installed on the evaporator that assists in the removal of the frost and ice. Although the heater helps speed up the frost and ice removal, it can raise the temperatures within the freezer.